

PROFESSIONAL SAAS USER GUIDE



Every Component. Every Movement. Every Time.

The Complete User Manual

Every module, workflow, and business invariant behind ImplaTrax's multi-tenant clinical inventory platform — from first login to production Supabase architecture.



Table of Contents

01 Welcome to ImplaTrax

Core philosophy, multi-tenant workspace isolation, and the append-only traceability principle.

02 Getting Started & Authentication

60-second onboarding, password vs. PIN, biometrics, and multi-device switching.

03 The Core Dashboard

Real-time KPI cards, trend charts, and global keyboard search.

04 Catalog & Product Taxonomy

Fixed vs. data-driven fields, SKU/barcode/QR generation, and bulk CSV import.

05 Inventory & Traceability Engine

The append-only movement ledger, manual adjustments, batch/lot tracking, and the stock ≥ 0 invariant.

06 Patients, Doctors & Clinical Cases

Doctor lookups, multi-case patients, and how placing a component creates a Sale atomically.

07 Sales Workflow

The mandatory patient rule, white-label invoices, and sale immutability.

08 Laboratory & Loans Workflow

The strict lab-only domain rule, partial returns, and lost-item accounting.

09 Vendors & Purchase Orders

The PO lifecycle, and the mandatory photo-evidence rule for every partial receipt.

10 Workspace Settings, Team & Security

White-label branding, least-privilege roles, device sessions, and the audit log.

This is an interactive PDF — every chapter above, and every cross-reference throughout the guide, is a clickable jump-link. Use your PDF reader's search (Ctrl/Cmd+F) to find any field, workflow, or business rule by name; every word in this document is real, selectable text.

CHAPTER 01

Welcome to ImplaTrax

The philosophy, architecture, and non-negotiable guarantees behind every screen in this product — read this chapter first, everything else builds on it.



CHAPTER 1

Welcome to ImplaTrax

ImplaTrax is a professional, multi-tenant inventory and traceability platform purpose-built for dental implant practices, labs, and the vendors who supply them. It exists to answer one question with total confidence, at any moment, for any component: where is it, who touched it, and what happened to it — from the day it was received to the day it was placed in a patient.

1.1 The Core Philosophy

"Every Component. Every Movement. Every Time."

That line is not a slogan bolted onto a marketing page — it is the design constraint every module in this guide is built around. Every unit of stock ImplaTrax knows about is traceable to a specific product, at a specific quantity, moving for a specific, recorded reason. Nothing enters or leaves inventory silently. Nothing is ever just "corrected" without a note explaining why. If a number on a dashboard looks wrong, the full history of exactly how it got there is always one click away — never reconstructed, never guessed.

1.2 Multi-Tenant Workspace Isolation

ImplaTrax runs on a shared, multi-tenant Supabase PostgreSQL database — every clinic, lab-facing account, and vendor-facing workspace lives in the same physical tables, but no workspace can ever see another workspace's data. This isn't a convention the application layer promises to respect; it's enforced by PostgreSQL Row Level Security (RLS) policies attached directly to every table, so even a bug in the front-end code, or a request crafted by hand against the API, cannot cross a workspace boundary.

- **Every business table carries a `workspace_id` column** — products, purchase orders, sales, loans, patients, cases, inventory movements, and every other record are all scoped to exactly one workspace.
- A `auth_workspace_ids()` policy function resolves which workspace(s) the currently signed-in user belongs to, and every RLS policy filters reads and writes against that set — a query simply cannot return rows outside it, regardless of what the client asks for.
- A person can belong to more than one workspace (e.g. a consultant supporting multiple clinics) and switches between them from the workspace picker — each workspace still keeps its own isolated products, team, PIN, and history.

**ARCHITECTURE INVARIANT****Isolation is enforced at the database, not the UI**

Row Level Security policies live on the PostgreSQL tables themselves. A workspace boundary in ImplaTrax is not "the app happens to only show you your own data" — it is "the database physically will not return another workspace's rows to your session, full stop."

1.3 The Append-Only Traceability Principle

The single most important architectural decision in ImplaTrax is that history is never rewritten. Several tables exist purely to record what happened, in order, forever:

- **Inventory Movements** — every stock change (a purchase order receipt, a sale, a loan going out, a loan returning, a manual adjustment, a lost item) writes one immutable row capturing the exact quantity before and after the change, who performed it, and why. Rows are never edited or deleted; a correction is a new movement, not a rewritten old one.
- **Purchase Order Events, Loan Events, Case Events** — every status transition on a purchase order, loan, or clinical case appends one entry to that record's own audit trail. A purchase order's history always shows "Created, then Submitted, then Partially Received, then Received" in that exact order, because each line was written once and never touched again.
- **The Workspace Audit Log** — security-relevant events (logins, PIN changes, member role changes, workspace renames) are appended the same way, giving every workspace a complete, tamper-evident compliance trail.

Practically, this means two things for anyone using the product: first, you can always answer "why does this product show 14 units on hand?" by reading its movement history top to bottom — the current quantity is nothing more than the running sum of every movement that ever touched it. Second, mistakes are corrected forward, never erased — voiding a sale or cancelling a purchase order adds new, clearly labeled history rather than pretending the original action never happened.

1.4 Who This Guide Is For

This manual documents ImplaTrax exactly as it behaves in production against the multi-tenant Supabase backend — every field, every business rule, every workflow described here is the real, enforced behavior of the system, not a simplified summary. It is written for three audiences at once: the clinical/front-desk staff who use ImplaTrax every day to receive stock, record sales, and manage loans; the workspace owners and admins responsible for team access and settings; and the engineers who need a precise reference for the invariants the system guarantees. Chapters 2 through 10 can be read independently once you've read this one — each documents a single functional area end to end, field by field, rule by rule.

CHAPTER 02

Getting Started & Authentication

From a clean-slate workspace to a fast, secure daily unlock — how ImplaTrax onboards a new team and keeps the operatory moving.




CHAPTER 2

Getting Started & Authentication

ImplaTrax separates two very different login moments: the one-time act of standing up a brand-new workspace, and the dozens of times a day a busy clinical team needs to get back into the app in under two seconds. Both are covered in this chapter.


2.1 60-Second Workspace Onboarding

A brand-new ImplaTrax workspace starts genuinely empty — there is no pre-installed product catalog, no sample vendors, and no placeholder patients to delete before you can trust what you see. Onboarding exists to get one real admin account and one real workspace name created, nothing more.

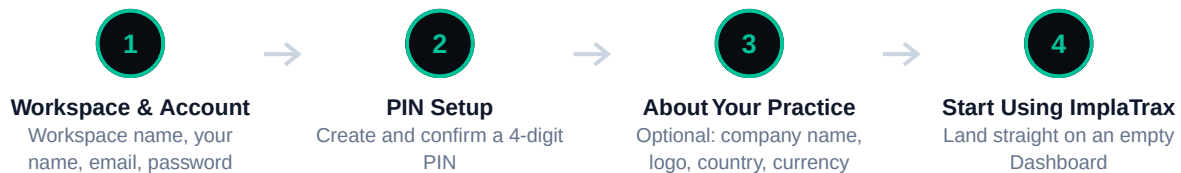

Every field on the "About Your Practice" step is explicitly optional — company name, workspace logo, country, and currency can all be filled in later from Settings without blocking the first login. The only two things onboarding requires are a workspace name and a working account (name, email, password) for the person setting it up, who becomes that workspace's **Owner**.

Onboarding — Account step

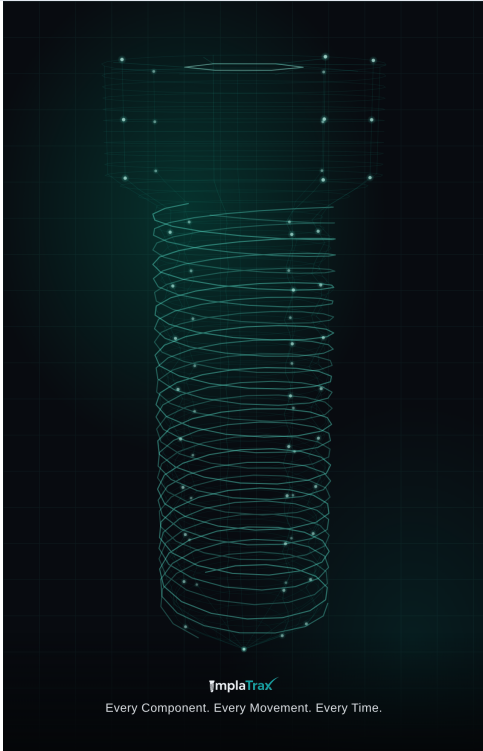



Every Component. Every Movement. Every Time.

Set Up Your Workspace

Tell us a little about your practice.

Workspace Name

Your Name

Email

Password

At least 8 characters. This is your account password — separate from the 4-digit PIN you'll set up next.

Continue

The account step collects the workspace name and the owner's credentials in a single, focused screen.

- 1 Workspace name — shown throughout the app and on every exported document.
- 2 Email + password — the account-level credential used to log in on a new device.
- 3 A 4-digit device PIN is created on the next step, for fast local unlock.

2.2 Two Ways to Sign In: Password vs. Device PIN

METHOD	BEST FOR	HOW IT WORKS
Email + Password	New device / new team member	The account-level credential. Used the first time a person signs in on a given device, or when a device has no PIN set for them yet.
4-Digit PIN	The daily clinical-operator unlock	A fast, device-scoped code — set once per device per person. A PIN created on one device is never valid on another; a new device always starts with a password sign-in.

This split exists because a clinical front desk cannot afford to type a full password every time the app locks between patients. The PIN is deliberately short and fast, but it is also deliberately device-scoped — stealing a PIN is useless anywhere except the one physical device it belongs to.

**SECURITY DESIGN****A PIN never travels between devices**

Each team member's PIN is tied to a specific device ID. Signing in on a new tablet or workstation always starts with the account password, followed by creating a brand-new PIN for that device — an old device's PIN cannot be guessed or reused elsewhere.

2.3 Biometric Authentication

Where the device supports it, ImplaTrax offers Touch ID, Face ID, or Windows Hello as a faster stand-in for the PIN, backed by the browser's WebAuthn platform authenticator API. Biometric unlock is opt-in per device and always has the PIN pad as a fallback — losing access to Face ID (a different person picking up the tablet, for instance) never locks anyone out of typing their PIN instead.

2.4 Multi-Device & Multi-Workspace Switching

- 1 A team member can be added to a workspace by an Owner, Super Admin, or Admin, who issues them an initial PIN directly — no email invitation loop required to get someone working on a shared device.
 - 2 A single email account can belong to more than one workspace (for example, a lab technician who supports two clinics) — the account menu's workspace switcher moves between them without signing out.
 - 3 Locking the app (manually, or automatically after inactivity) always returns to the PIN pad for whoever is signed in — it never signs the whole device out to the workspace picker unless "Log Out" or "Switch User" is chosen explicitly.
- **Auto-lock** — configurable minutes of inactivity before the app locks itself (0 disables it).
 - **Session timeout** — an absolute cap on how long a session stays valid since the last unlock, regardless of activity.
 - **Failed-attempt lockdown** — a maximum number of incorrect PIN attempts triggers a temporary, timed lockdown before another attempt is allowed.

All three of these are workspace-wide security preferences, configurable by a workspace manager from Settings, under Security — see Chapter 10.

CHAPTER 03

The Core Dashboard

One screen, six real-time signals — what changed, what's low, and what's still open, the moment you log in.



CHAPTER 3

The Core Dashboard

The Dashboard is deliberately the first thing every workspace member sees after unlocking ImplaTrax — a live snapshot of the numbers that matter, computed directly from the same append-only data every other chapter of this guide describes, never a cached or delayed summary.

3.1 Real-Time KPI Cards

Six cards sit across the top of the dashboard, each reading directly from current workspace state:

Inventory Value

Total on-hand stock valued at unit cost across every product — the true replacement value of what's on the shelf right now.

Low Stock Items

A live count of active products at or below their configured reorder threshold — the same set that drives the sidebar's low-stock alerts.

Open Loans

Loans currently issued to labs that are not yet fully returned or closed — your outstanding B2B exposure at a glance.

Pending Purchase Orders

Draft, submitted, confirmed, and partially-received orders — everything still short of a full receipt.

Cases This Month


Clinical cases created in the current calendar month, a fast read on chairside volume.

Revenue (30d)

Total value of non-voided sales recorded in the trailing 30 days.

3.2 Trend Charts

- **Stock Movements** — inbound vs. outbound quantity over the last 14 days, drawn straight from the inventory movement ledger described in Chapter 5.
- **Inventory Value by Manufacturer** — current stock valuation broken out across your product catalog's manufacturer field, useful for spotting where capital is concentrated.

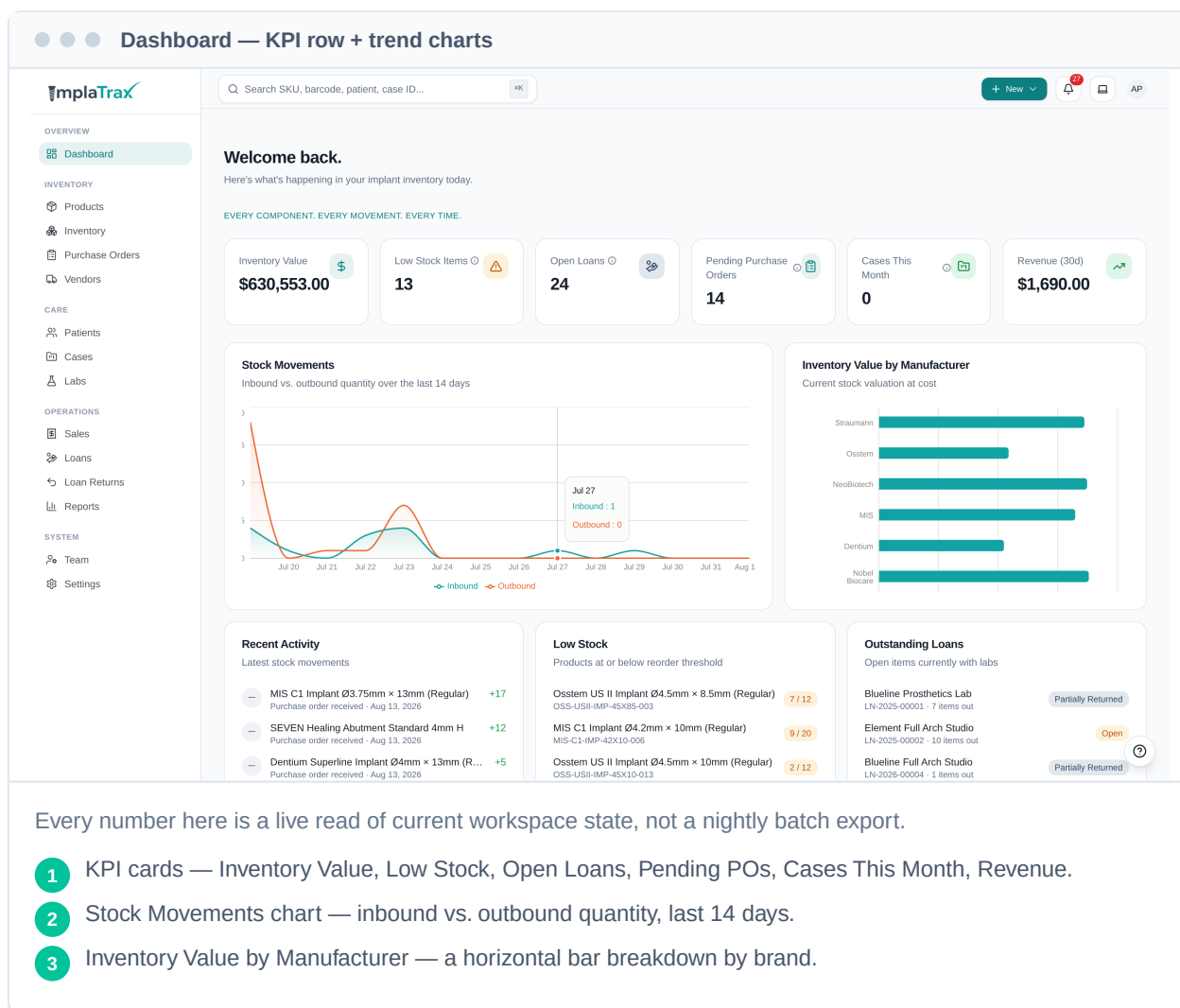

3.3 Global Keyboard Search

Press **Ctrl+K** (or **Cmd+K** on macOS) from anywhere in ImplaTrax to open Global Search — a single command palette that searches across product names, SKUs, barcodes, patients, and clinical case IDs simultaneously, and jumps straight to the matching record. It's the fastest way to answer "do we have this part" or "which case is this component tied to" without leaving whatever screen you're currently on.



TIP

Barcode scanners work here too

Any USB or Bluetooth barcode scanner configured as a keyboard-emulation device can type a scanned barcode straight into Global Search — scan a physical component and ImplaTrax opens its product record instantly, no dedicated scanning mode required.

3.4 Activity, Low Stock & Loan Widgets

- Recent Activity** — the newest inventory movements across the whole workspace, most recent first.

- **Low Stock** — every product at or below its reorder threshold, ranked so the most depleted items surface first.
- **Outstanding Loans** — open and partially-returned loans, so a lab relationship never quietly goes stale.

CHAPTER 04

Catalog & Product Taxonomy

How ImplaTrax classifies every implant, abutment, and consumable — and what changes per product versus platform-wide.



CHAPTER 4

Catalog & Product Taxonomy

Every product ImplaTrax tracks is described the same way, whether it was typed in one at a time or arrived in a 500-row import — a consistent shape is what makes search, barcode lookup, and reporting reliable across an entire catalog.

4.1 Fixed Taxonomy vs. Free-Text Fields

Two of a product's classification fields — **Manufacturer** and **Category** — come from a fixed, platform-wide list (six real-world implant manufacturers; ten component categories from Implant Fixture through Prosthetic Screw). This is deliberate: these are facts about the dental implant industry, not something one workspace should be able to define differently from another, so every workspace's reports and filters mean the same thing. **System** and **Platform**, by contrast, are free-text fields captured per product (e.g. "BLX", "TS III", "TiBase" for system; "NC", "RC", "WP" for platform) — genuinely data-driven, since implant systems and platform codes vary by manufacturer far more than can be usefully enumerated in a fixed list.

FIELD	TYPE	DESCRIPTION
manufacturer	fixed list	One of six real implant brands — the same list every workspace shares.
category	fixed list	One of ten component categories (Implant Fixture, Healing Abutment, Cover Screw, Bone Graft Material, etc).
system	free text	The manufacturer's implant line, e.g. "BLX", "TS III" — captured per product.
platform	free text, optional	Connection platform code, e.g. "NC" / "RC" / "WP" — varies by manufacturer and product.
diameterMm / lengthMm	number, optional	Physical dimensions, shown wherever a product is listed so the right size is never guessed.

4.2 Identity: SKU, Barcode & QR

Three identifiers are generated automatically the moment a product is created and never need to be typed by hand:

- **SKU** — composed from the manufacturer, system, and a sequence number (e.g. STR-BLX-014), unique within the workspace.
- **Barcode** — a CODE128-compatible value, rendered as a scannable barcode on every printed document and product label.


- **QR Payload** — encodes the product's SKU and internal ID, giving a phone camera a fast path straight to the correct product record.

4.3 Pricing & Visibility

FIELD	DESCRIPTION
unitCost	What the workspace pays the vendor — drives Inventory Value on the Dashboard.
unitPrice	What the workspace charges the patient — used as the default line price on a new Sale.
priceVisible	A per-product toggle for whether pricing shows to staff by default; the workspace-wide "Price Visibility" setting in Clinic Settings supplies the default for new products.
lowStockThreshold	The quantity at or below which this product is flagged as low stock — feeds the Dashboard, the sidebar alert, and Reports.
batchTracked	Whether this specific product participates in lot/batch capture — only meaningful when the workspace-wide Batch/Lot Tracking setting is on (see Chapter 5).

4.4 Bulk CSV Import

Migrating from a spreadsheet, or seeding a large catalog in one pass, goes through the CSV importer rather than one-at-a-time product creation.

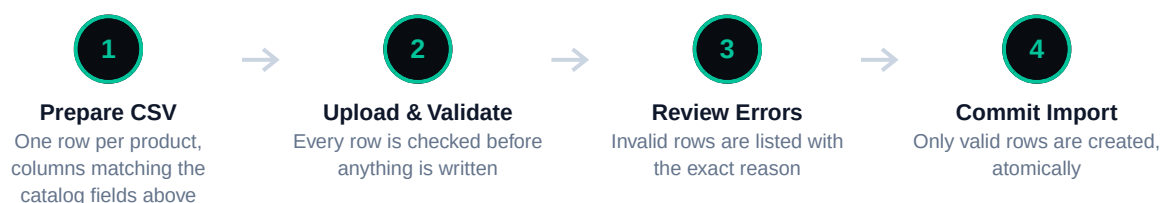

DATA INTEGRITY

A bad row never blocks a good one — but validation runs first

The importer validates every row up front and shows exactly which rows would fail and why, before anything is written to the catalog. You choose to fix the file and re-upload, or proceed with only the valid rows — there is no partial-write state where some rows silently succeed while others silently fail.

- 1 Manufacturer and Category values in the file must match Implastrax's fixed lists exactly (case-insensitive) — anything else is flagged as an error rather than silently dropped or guessed.

- 2 Newly imported products immediately get their SKU, barcode, and QR payload generated the same way a manually created product does.
- 3 If any imported product starts with a non-zero quantity, that quantity is recorded as an "Initial stock on bulk import" inbound movement — the append-only ledger never has a product appear with stock and no explanation for it.

CHAPTER 05

Inventory & Traceability Engine

The append-only movement ledger that makes every unit of stock explainable, and the invariant that keeps it honest: stock can never go negative.



CHAPTER 5

Inventory & Traceability Engine

Every quantity ImplaTrax shows you — a product's stock count, a workspace's total inventory value — is a derived number, never a manually maintained one. This chapter explains what actually produces it.

5.1 The Append-Only Movement Ledger

Every single event that changes how much of a product is on hand — a purchase order receipt, a sale, a loan going out, a loan returning, a manual correction, a component reported lost — writes exactly one row to the `inventory_movements` table. A product's current quantity is nothing more than the running sum of its movements; it is never stored or edited as an independent number that could drift out of sync with its own history.

FIELD	TYPE	DESCRIPTION
type	<code>enum</code>	inbound · outbound · adjustment · loan-out · loan-return · sale · lost
quantity	<code>signed integer</code>	Positive = stock increase, negative = decrease — the exact delta this one event applied.
quantityBefore / quantityAfter	<code>integer</code>	A full snapshot of the product's stock immediately before and after this movement — self-contained, never recomputed from other rows.
reason	<code>text</code>	Always present. System-generated for automated events ("Purchase order received"), required and operator-typed for manual adjustments.
reference	<code>text, optional</code>	The human-readable business document this movement belongs to — a PO number, Sale number, or Loan number.
performedBy	<code>user id</code>	Exactly who triggered this movement — attributed at write time, not reconstructed later.



ARCHITECTURE INVARIANT

Movements are never edited or deleted

Correcting a mistake always means writing a new, clearly-reasoned movement — never editing or removing an old one. A product's full movement history is a permanent, chronological, tamper-evident record of everything that ever happened to it.

5.2 Manual Stock Adjustments

Physical counts sometimes disagree with what ImplaTrax expects — damaged stock, a miscount, an item found in the wrong bin. The Adjust Stock action exists for exactly this, and it enforces one rule strictly:

**MANDATORY FIELD****A reason is required for every manual adjustment**

ImplaTrax rejects a manual stock adjustment with an empty or whitespace-only reason — this is enforced in the data layer itself, not only the form, so there is no code path that can silently correct a quantity without an explanation attached to it forever.

An adjustment can move quantity in either direction, but the resulting on-hand quantity is always clamped at zero — see the invariant below.

5.3 Batch/Lot Tracking & Expiry

Batch/Lot Tracking is a single workspace-wide switch in Clinic Settings — when it's on, every purchase order receipt requires a lot number for the line being received, regardless of whether that specific product is individually flagged as batch-tracked. Traceability, once it's turned on, begins at the moment stock enters the workspace and is never optional on a per-receipt basis.

- Each lot number received creates its own **Product Batch** record, carrying the quantity received, an optional expiry date, and the purchase order it came from.
- Receiving the **same** lot number across two separate deliveries creates two batch records, not one merged one — every receipt keeps its own point-in-time trail.
- The lot a component carries stays attached through Sales and Loans, so a specific implant placed in a specific patient can always be traced back to the exact delivery it arrived on.

5.4 Barcode Lookup

A product's barcode is a first-class, indexed field, not just something printed on a label — typing or scanning a barcode value into Global Search (Chapter 3) jumps straight to that product's record, which is the fastest way to confirm identity and current stock before it's committed to a sale or loan.

● ● ●

Product detail — movement history

Nobel Biocare NobelParallel CC

Implant Ø5mm x 10mm (NP)

© NOB-NBIE-IMP-SX10-009

Active

Normal

Nobel Biocare

Implant Future

On hand

72

Available

72

Reserved

0

Reorder level

13

Vendor

Nobel Biocare Distribution

Unit cost

\$111.00

Unit price

\$183.00

Titanium root-form dental implant fixture by Nobel Biocare, engineered for primary stability across a range of bone densities.

Identifiers

Barcode

899310000000

QR Code

Manual stock adjustment

-

1

+

Reason (required)

e.g. Cycle count adjustment

This note becomes part of the permanent inventory history.

Remove stock

Add stock

History

Purchase History

DentalSource Wholesale

Dec 11, 2025, 3:59 PM SL-2025-12041559

+8

Sales History

Lisa Campbell

Feb 26, 2025, 3:59 PM SL-2025-00028

-2

Loan History

Richard Flores - Dr. Alan Whitfield

Jul 1, 2025, 3:59 PM SL-2025-00018

-1

Loan History

Artisan Dental Arts

Jul 19, 2025, 3:59 PM LN-2025-00013

-6

Loan History

Artisan Dental Arts

Jul 12, 2025, 3:59 PM LN-2025-00013

+8

Loan History

Menden Implant Lab

Mar 31, 2025, 3:59 PM LN-2025-00014

-3

Loan History

Menden Implant Lab

Mar 2, 2025, 3:59 PM LN-2025-00014

+3

Loan History

Horizon Dental Studio

Dec 18, 2025, 3:59 PM LN-2025-00012

-6

Adjustments

No manual adjustments recorded yet.

Every row on this timeline is an immutable inventory_movements entry — the running total at the top is simply their sum.

- Current quantity on hand, low-stock threshold, and batch-tracking status.
- Movement ledger — type, signed quantity, before/after snapshot, reason, and who performed it.

5.5 The Stock ≥ 0 Invariant

DATABASE INVARIANT

Stock can never drop below zero

quantity_on_hand carries a database-level CHECK (quantity_on_hand ≥ 0) constraint. No sale, loan, or adjustment — however it's requested — can push a product's stock negative; the operation is rejected (or, for a downward adjustment, the resulting quantity is clamped at zero) before it can ever leave inventory in an impossible state. Available-stock checks in the UI exist to give a friendly warning before that point, not as the only line of defense.


ImplaTrax Professional SaaS User Guide — Version 1.0

Page 21 of 39

CHAPTER 06

Patients, Doctors & Clinical Cases

Where inventory meets the clinic — how a case links a patient, a doctor, and the components actually placed.

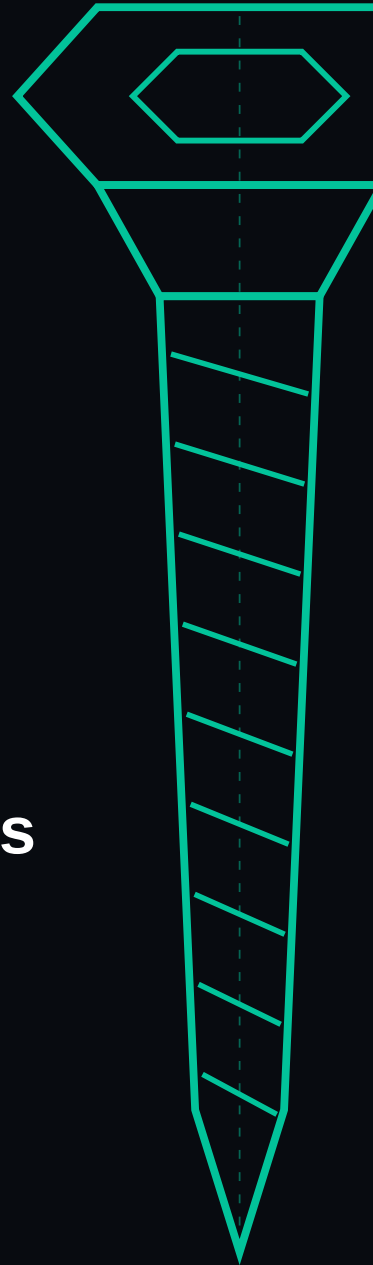

CHAPTER 6

Patients, Doctors & Clinical Cases

A Case is where ImplaTrax's inventory engine and clinical record-keeping meet — it is the record that proves which specific components ended up in which specific patient, under which doctor's care.

6.1 Doctors

Doctors are a real, workspace-level lookup table — not a free-typed name repeated on every patient and case. A doctor is typically created inline, the first time their name is needed on a patient or case form, and from then on is reused everywhere via a searchable picker. This keeps "Dr. Alan Whitfield" spelled identically across every record he's associated with, which matters the moment you need to report on a single doctor's case volume or component usage.

- Doctors are never hard-deleted — a doctor no longer active at the practice is archived instead, disappearing from the picker for new selections while every existing patient/case record that already references them stays completely intact.


6.2 Patients

FIELD	TYPE	DESCRIPTION
patientCode	auto-generated	A short, human-readable identifier (e.g. PT-01042), assigned the moment the patient is created.
firstName / lastName / dob / sex	text / date / enum	Core identity fields for matching a patient record to the person in the chair.
phone / email	text, optional	Contact details, shown on the patient profile and searchable from Global Search.
primaryDoctor	text	The patient's usual treating doctor — a display string sourced from the Doctors lookup, distinct from any doctor recorded on an individual case.

A patient's profile aggregates every case ever opened for them, so a coordinator can see a person's complete implant history — not just the case currently in progress — from one screen.

6.3 Clinical Cases

A single patient can have multiple cases running in parallel — a full-arch case and a separate single-tooth case, for instance — each tracked completely independently: its own status, its own doctor, its own list of components placed.

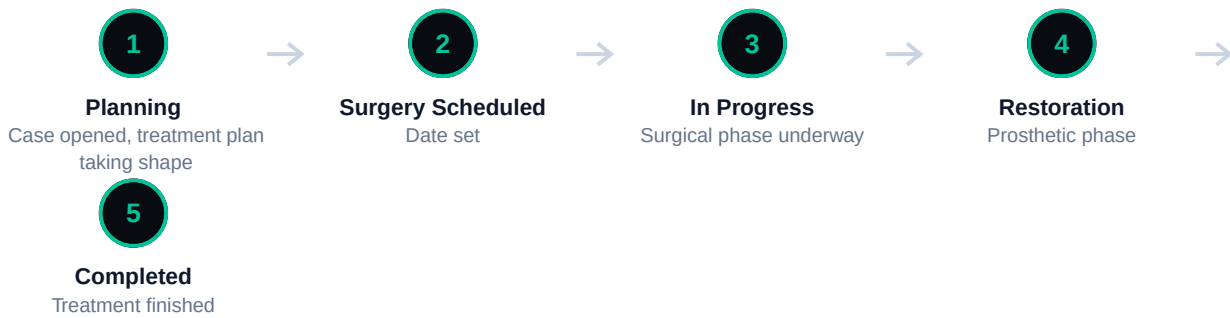

Every case also carries an append-only history — status advances, and each implant added — recorded the same way a purchase order's or loan's history is (Chapter 1). A case can also be cancelled at any stage instead of progressing to completion; both Completed and Cancelled are terminal states the guard logic will not allow any further status transition out of.

6.4 Adding a Component to a Case Is a Sale

Recording that an implant, abutment, or other component was placed in a patient during a case does not just add a line to a notes field — it atomically deducts the component from stock, writes an **outbound** inventory movement, and creates a real Sale record linked to that patient and case, all as one indivisible action.



DATA INTEGRITY

A case's implant list can never drift out of sync with stock or sales

If the underlying Sale would fail for any reason — insufficient stock, for instance — the whole action is rejected before anything is written: the case's implant list, the product's stock, and the sales ledger either all update together or none of them do.

Every implant recorded on a case carries the tooth number it was placed at (FDI notation, e.g. "36") and, when Batch/Lot Tracking is on, the exact lot it came from — giving a case a complete, component-level traceability record from delivery to placement.

CHAPTER 07

Sales Workflow

Clinical consumption, the mandatory patient rule, and the white-label documents every sale can generate.



CHAPTER 7

Sales Workflow

A Sale in ImplaTrax means one thing, unconditionally: a component was permanently placed in — or delivered to — a real patient. It is the system's clinical-consumption record, and it is built around one non-negotiable rule.

7.1 The Mandatory Patient Rule




CRITICAL BUSINESS RULE

A sale cannot exist without a patient

ImplaTrax has no concept of an anonymous or walk-in sale. Every Sale record requires a valid `patientId` — the action is rejected outright, before any stock is touched, if one isn't supplied. This is enforced in three independent places: the Record Sale form won't submit without a patient selected, the data layer rejects the action even if called directly, and the database schema itself makes `sales.patient_id` a not-null, foreign-key-constrained column. There is no path — UI, API, or direct database write — that produces a sale disconnected from a real patient.

This exists because a Sale is ImplaTrax's medical traceability record — the whole point of the platform is being able to answer, with certainty, exactly which patient received a specific implant lot. A sale with no patient would be a permanent hole in that record.

7.2 Recording a Sale, Step by Step

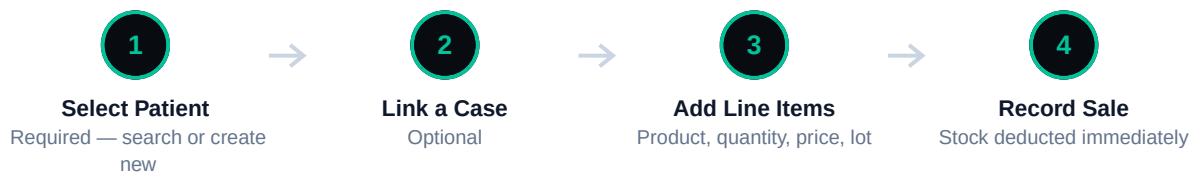

- 1 Choose the patient the sale is for — an existing patient found via search, or a brand-new one created inline without losing any line items already entered.
- 2 Optionally link the sale to one of that patient's existing clinical cases — sales linked to a case are labeled "Used in patient case" in the movement history; unlinked sales are labeled "Direct sale".
- 3 Add one or more product lines, each with its own quantity and unit price (defaulting to the product's configured price, editable per line) — a batch/lot field appears automatically when Batch/Lot Tracking is on for a batch-tracked product.
- 4 Every line's requested quantity is checked against current stock, combined across lines for the same product, before the sale can be submitted — you'll see exactly how many units are actually available if a line exceeds it.

- Submitting the sale deducts stock for every line, writes one 'sale' movement per line, and generates the sale's number (e.g. SL-2026-00042) — all as a single atomic action.

7.3 White-Label Document Generation

Every sale can produce a printable Sales Invoice and Delivery Challan, and both documents are built around a deliberate brand hierarchy: your workspace's own identity is the hero, Implatrax is a quiet footnote.

Sales Invoice — white-label header

INVOICE NO.

SL-2026-00001

INVOICE DATE

Mar 7, 2026

PATIENT

Matthew Jackson

CASE

IDC-2025-00008

SOLD BY

Unknown

PRODUCT	QTY	UNIT PRICE	LINE TOTAL
Osstem TS III Implant Ø5mm × 10mm (Wide) OSS-TSIII-IMP-5X10-070	1	\$261.00	\$261.00
Straumann TLX Implant Ø4.5mm × 10mm (TL) STR-TLX-IMP-45X10-001	1	\$271.00	\$271.00

AMOUNT DUE

\$532.00

This is a computer-generated document from Implatrax.
Generated Aug 1, 2026, 3:59 PM

POWERED BY 

1

Workspace logo + clinic name/address — the primary, large-format brand.

2

Invoice number, patient, case reference, and line items.

3

"Powered by Implatrax" — small footer attribution, ~25% the size of the workspace logo.

- If no workspace logo has been uploaded yet, the header falls back gracefully to a styled clinic name — the layout never shows an empty gap where a logo would go.
- A workspace logo uploaded or changed in Settings (Chapter 10) applies to every document generated from that point forward, immediately — no separate publish step.

7.4 Sale Immutability & Reversal

A completed sale's line items, quantities, and prices are never edited after the fact — a sale is a permanent clinical record, not a draft. If a sale genuinely needs to be undone (an item was recorded in error, a patient declined treatment after checkout), the reversal path is **Void**, not edit or delete.

**APPEND-ONLY DESIGN****Voiding restores stock without erasing the record**

Voiding a sale requires a typed reason, restores every line's quantity back to stock, and marks the sale as voided with a timestamp and that reason — the original sale record stays fully visible for audit purposes, permanently excluded from revenue totals but never deleted. History is corrected forward, never rewritten (see Chapter 1).

CHAPTER 08

Laboratory & Loans Workflow

Strictly B2B: components on loan move only between the workspace and dental laboratories — never patients, never vendors.



CHAPTER 8

Laboratory & Loans Workflow


A Loan tracks components that leave the workspace temporarily — surgical kits, trial components, consignment stock sent to a lab for a specific case — with an expectation that most of it comes back.

8.1 Strict Domain Rule: Labs Only

**CRITICAL BUSINESS RULE****A loan can only be issued to a lab**

Every Loan requires a valid `labId` referencing a real lab record in the workspace — the action is rejected if the ID doesn't resolve to one. Loans exist to model a specific, recurring B2B relationship (a clinic and the dental laboratories it works with); they are not a general-purpose "lend to anyone" mechanism. A patient or vendor is never a valid loan counterparty — component movement to a patient is always a Sale (Chapter 7), and movement from a vendor is always a Purchase Order receipt (Chapter 9).

8.2 Issuing a Loan

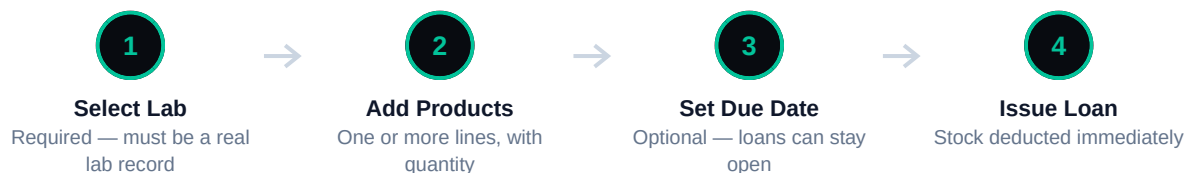

Issuing a loan deducts every line's quantity from stock right away and writes a **loan-out** movement per line, referencing the loan's number (e.g. LN-2026-00012). Like sales, the combined requested quantity for a product across multiple lines is checked against real stock before the loan can be issued — and when Batch/Lot Tracking is on, a batch-tracked product requires a lot number captured at issuance, since traceability begins the moment stock leaves the building, not only when it returns.

8.3 Tracking Open Loans & Processing Returns

LOAN STATUS	MEANING
Open	Issued, nothing returned or reported lost yet.
Partially Returned	Some, but not all, of the loaned quantity has been accounted for.
Closed	Every unit loaned has been either returned or reported lost — nothing left outstanding.

A return is processed per line: enter how many units came back, and separately, how many are being reported lost. Returned quantity restores directly to stock and writes a **loan-return** movement; lost quantity does **not** restore stock — it was already deducted at issuance — but is recorded with its own **lost** movement documenting the disposition.

**MANDATORY FIELD****A reason is required for every lost component**

Reporting any quantity as lost on a return is rejected without a typed reason — every unit that doesn't come back has a documented explanation attached to it permanently, the same discipline applied everywhere else stock is corrected outside the normal flow.

- 1 Partial returns are fully supported — a loan can be returned against multiple times over its life, each partial return appending its own entry to the loan's history.
- 2 A loan automatically moves to Closed the moment returned + lost quantity accounts for every unit originally loaned, across every line.
- 3 A closed loan cannot be returned against again — attempting to do so is rejected by the same guard logic that governs every other terminal-state record in ImplaTrax.

8.4 White-Label Loan Documents

Both the initial Loan Challan (what was issued) and the Return Delivery Slip (what came back) follow the same white-label hierarchy as Sales documents in Chapter 7 — the workspace's own logo and identity lead the header, with ImplaTrax reduced to a small "Powered by" footer mark.

- Optional photo attachments can be added at loan issuance and again at return — useful for recording component condition or a signed physical delivery slip — and are linked to the specific inventory movement(s) that action created.

CHAPTER 09

Vendors & Purchase Orders

From draft to received stock — including the mandatory photo-evidence rule that protects every partial delivery.



CHAPTER 9


Vendors & Purchase Orders

Stock enters ImplaTrax through exactly one door: a Purchase Order being received. This chapter covers vendor profiles, the full PO lifecycle, and the mandatory-photo rule that governs every partial delivery.

9.1 Vendor Profiles

A vendor record holds contact details, country, and — critically — the set of manufacturers that vendor supplies. That manufacturer association is what lets a purchase order line suggest the right vendor for a given product automatically, and lets Reports break spend down by supplier relationship rather than just by raw product.

9.2 The Purchase Order Lifecycle

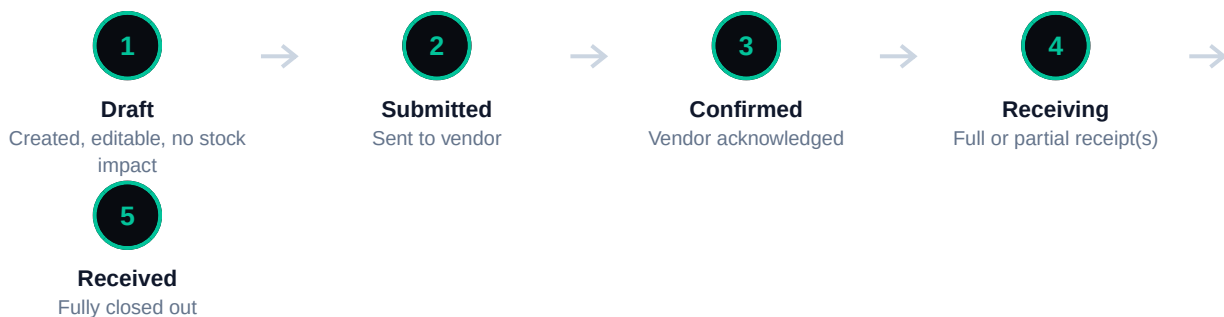

ARCHITECTURE INVARIANT

A purchase order never touches inventory by itself

Creating, submitting, or confirming a purchase order only ever changes that order's own status and append-only history — never a product's stock. Inventory changes exactly once a receipt is recorded against it, and not a moment before. A draft PO for 500 units has zero effect on what the Dashboard reports as on-hand.

A purchase order can also be **cancelled** — but only before it has been fully received; once every line is received, the order is closed out and can no longer be cancelled, since doing so would contradict inventory that has already physically arrived.

9.3 Receiving Stock: Full vs. Partial

Receiving is where a purchase order actually affects inventory. For each line on the order, you enter how many units arrived — up to, but never exceeding, what's still outstanding for that line. A receipt is treated as **partial** the moment any line's received quantity comes in short of what was outstanding for it — whether

that line was under-received this time, or skipped entirely while other lines on the same order were received in full. A receipt only counts as **full** when every line on the order is completely caught up.

FIELD	DESCRIPTION
Full Receipt	Every line's received quantity now equals what was ordered. Order status moves to Received; photo evidence is optional.
Partial Receipt	At least one line remains short of its ordered quantity after this receipt. Order status moves to Partially Received; photo evidence is required — see 9.4.

9.4 Mandatory Photo Evidence for Partial Receipts



CRITICAL BUSINESS RULE

A partial receipt cannot be submitted without a photo

Whenever a receipt is partial, ImplaTrax requires at least one photo of the delivery slip or the physical package before the receipt can be confirmed. The Confirm Receipt button stays disabled and a clear inline error is shown until a photo is attached; the same rule is enforced independently by the data layer — a partial receipt submitted with zero photos throws and nothing about the order or inventory is changed. A full receipt is never blocked this way; photos are optional once every line is fully caught up.

The reasoning is straightforward: a partial delivery is exactly the moment a discrepancy — short-shipped units, a damaged carton, a substitution — is most likely, and hardest to reconstruct later from memory. Requiring photographic evidence at the moment of receipt turns "the vendor says they sent 50, we only got 40" from a dispute into a documented fact.

● ● ● Receive Purchase Order — partial receipt blocked

×

Receive 202602181559

Enter the quantity received for each line item. Partial receiving is supported.

IS-III Active Healing Abutment Regular 6mm H

Ordered 55 · Received 0 · Remaining 55

27

CopiOs Resorbable Membrane 20×30mm

Ordered 21 · Received 0 · Remaining 21

21

CopiOs Resorbable Membrane 20×30mm

Ordered 16 · Received 0 · Remaining 16

16

IS-III Active Impression Coping — Open Tray Reg...

Ordered 17 · Received 0 · Remaining 17

17

CopiOs Bone Graft 1.0cc

Ordered 56 · Received 0 · Remaining 56

56

Shipment Photos (Required for Partial Receipt) *

Add

Please attach at least one photo of the delivery slip or package to document this partial receipt.

Cancel

Confirm Receipt

Entering fewer units than remain outstanding on any line flips the photo field to required; Confirm Receipt stays disabled until at least one photo is attached.

- 1 Per-line quantity entry, capped at what remains outstanding for that line.
- 2 "Shipment Photos (Required for Partial Receipt)" — label and inline error change dynamically.
- 3 Confirm Receipt — disabled until the photo requirement is satisfied.

- 1 Photos attached to a receipt are linked to every inventory movement that receipt creates, and accumulate — append-only — on the purchase order itself across however many partial receipts it takes to fully close it out.
- 2 When Batch/Lot Tracking is on, a lot number is required for every line being received in the same action — the photo requirement and the lot requirement are independent checks, both must pass to submit.
- 3 Receiving a batch-tracked line with a lot number creates a Product Batch record and stamps that lot onto the resulting inventory movement, exactly as described in Chapter 5.

9.5 Automatic Restocking & Movement History

Confirming any receipt — full or partial — immediately increases the affected products' on-hand quantity and writes one **inbound** movement per line received, referencing the purchase order's number. There is no separate "apply to inventory" step; receiving and restocking are the same atomic action.

ImplaTrax Professional SaaS User Guide — Version 1.0

Page 35 of 39

CHAPTER 10

Workspace Settings, Team & Security

Brand identity, least-privilege team access, device sessions, and the audit log that ties it all together.



CHAPTER 10

Workspace Settings, Team & Security

Settings is where a workspace's identity, its team's access, and its security posture all live — every change here is scoped to your workspace alone, protected by the same Row Level Security boundary described in Chapter 1.

10.1 Workspace Name & White-Label Logo

The Workspace Name field is shown throughout the app and on every exported document. A dedicated Workspace Logo card lets any workspace manager upload, change, or remove the workspace's logo at any time — independent of whether one was set during onboarding.

● ● ● Settings / Workspace / Workspace Logo

 **Workspace Logo**

Appears as the primary brand on invoices, challans, and reports you export or print



Change Logo

✕ Remove

PNG or JPG, square works best. Optional — falls back to your workspace name if none is set.

Upload, Change, and Remove — the logo syncs live across the app the instant it changes, no save step or page reload required.

1 Live preview of the current logo, or a clean placeholder when none is set.

2 Upload / Change — replaces the logo instantly.

3 Remove — reverts every surface to the styled workspace-name fallback.

- A newly uploaded logo appears immediately on every printed Sale, Loan, and Purchase Order document generated afterward — see the white-label hierarchy described in Chapters 7–9.

- The workspace logo is intentionally scoped to Settings and printed documents only — the in-app navigation shell (sidebar and top bar) always shows the ImplaTrax mark alone, so the product's own identity stays unambiguous while you're using it.

10.2 Team Members & Least-Privilege Roles

ROLE	CAN MANAGE WORKSPACE/TEAM?
Owner	Full access — every workspace has exactly one, the account that completed onboarding.
Super Admin / Admin	Full workspace and team management — add/remove members, change roles, rename the workspace, upload the logo.
Manager	Day-to-day operational access without workspace-level administration rights.
Staff	Standard clinical/operational access.
Read Only	View access across the workspace, no create/edit/delete actions.

A team member is added by an Owner, Super Admin, or Admin, who issues their initial PIN directly — the fastest path to a new hire being productive on a shared clinic device the same day, without waiting on an email invitation round-trip.



SAFEGUARD

A workspace can never lose its last active Owner

ImplaTrax will not let the final active Owner be disabled, removed, or demoted while they're the only one left — a workspace can never accidentally end up with no one able to manage it.

10.3 Device Sessions & PIN Resets

- **Reset PIN** — a manager can force any team member to set a brand-new PIN on their next unlock, useful if a device is lost or a PIN is suspected compromised.
- **Disable / Reactivate** — revokes or restores a member's access without deleting their history — every record they ever touched stays attributed to them.
- **Sessions & Devices** — a per-device record of this browser install (a stable device ID), used to scope PINs and show recent activity.

Three workspace-wide security preferences round out access control, all configurable from Settings, under Security:

FIELD	DESCRIPTION

Auto-Lock	Minutes of inactivity before the app locks itself and returns to the PIN pad. Zero disables it.
Session Timeout	An absolute cap, in minutes, on how long a session stays valid since the last unlock — regardless of activity.
Lockout	Maximum consecutive failed PIN attempts before a timed lockout, and how long that lockout lasts.

10.4 The Append-Only System Audit Log

Every security-relevant event — logins, lockouts, PIN changes and resets, biometric enrollment, member additions/removals/role changes, workspace renames — is appended to the workspace's Activity Log the same way every other history table in ImplaTrax works: written once, in order, never edited or removed.



COMPLIANCE

A complete, tamper-evident trail for every workspace

Because the audit log is append-only and scoped by the same Row Level Security boundary as every other table, it gives each workspace an accurate, complete record of who did what and when — without any workspace being able to see, or worry about, another workspace's activity.